


Melker Cuffed Emergency Cricothyrotomy Catheter Sets

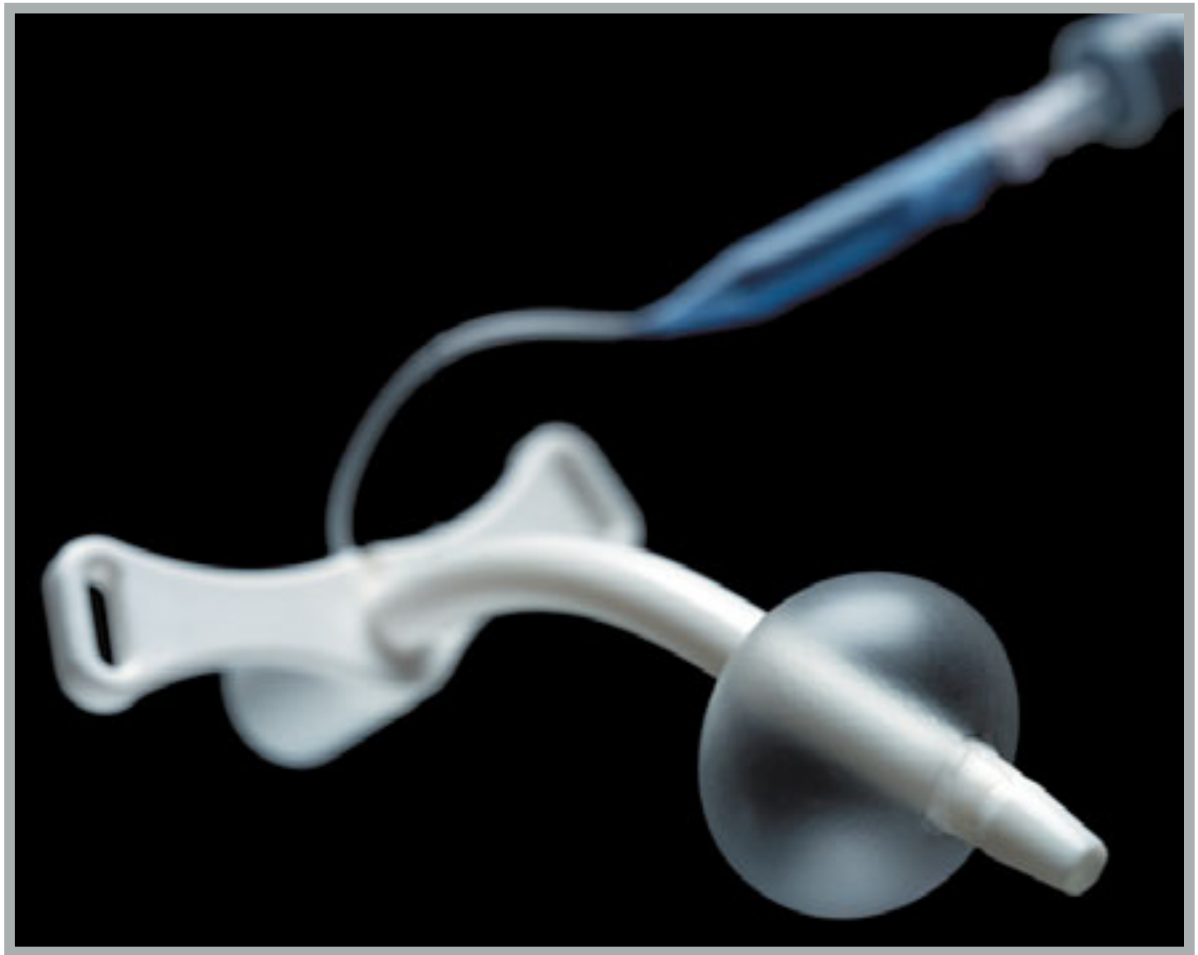


COOK[®]

An emergency airway catheter that enables control and protection of the airway.

INFLATED CUFF

When inflated, the cuff creates an airway seal to prevent aspiration and control ventilation.

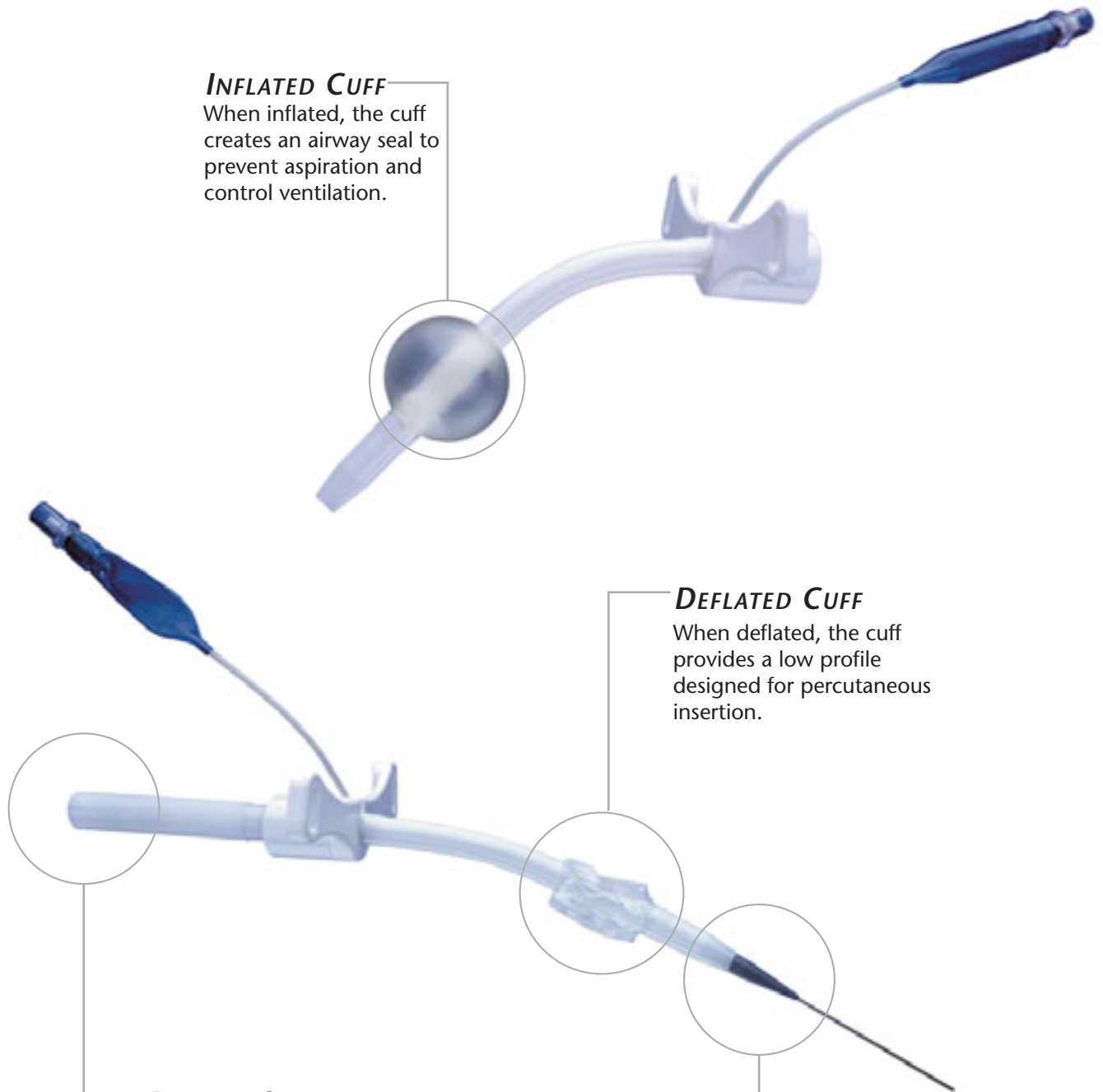

DEFLATED CUFF

When deflated, the cuff provides a low profile designed for percutaneous insertion.

PROVEN CONTROL


Sturdy, yet flexible enough for rapid insertion, the contoured handle helps prevent slippage.

SELDINGER TECHNIQUE

Engineered to be inserted via the Seldinger Technique. The use of this over the wire method is a recognized standard for percutaneous entry.

MELKER EMERGENCY CRICOTHYROTOMY CATHETER SETS


CUFFED

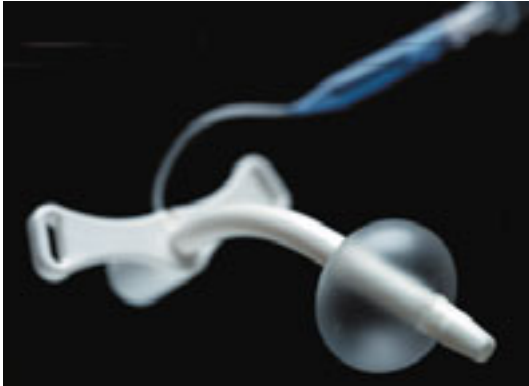

Used for emergency airway access when conventional endotracheal intubation or ventilation may not be possible.

- Airway access is achieved utilizing percutaneous entry (Seldinger) technique via the cricothyroid membrane.
- Subsequent dilation of the tract and tracheal entrance site permits passage of the emergency airway.
- Balloon cuffed catheter protects and controls the airway.
- Airway catheter has radiopaque stripe.

Supplied sterile in peel-open packages. Intended for one-time use.

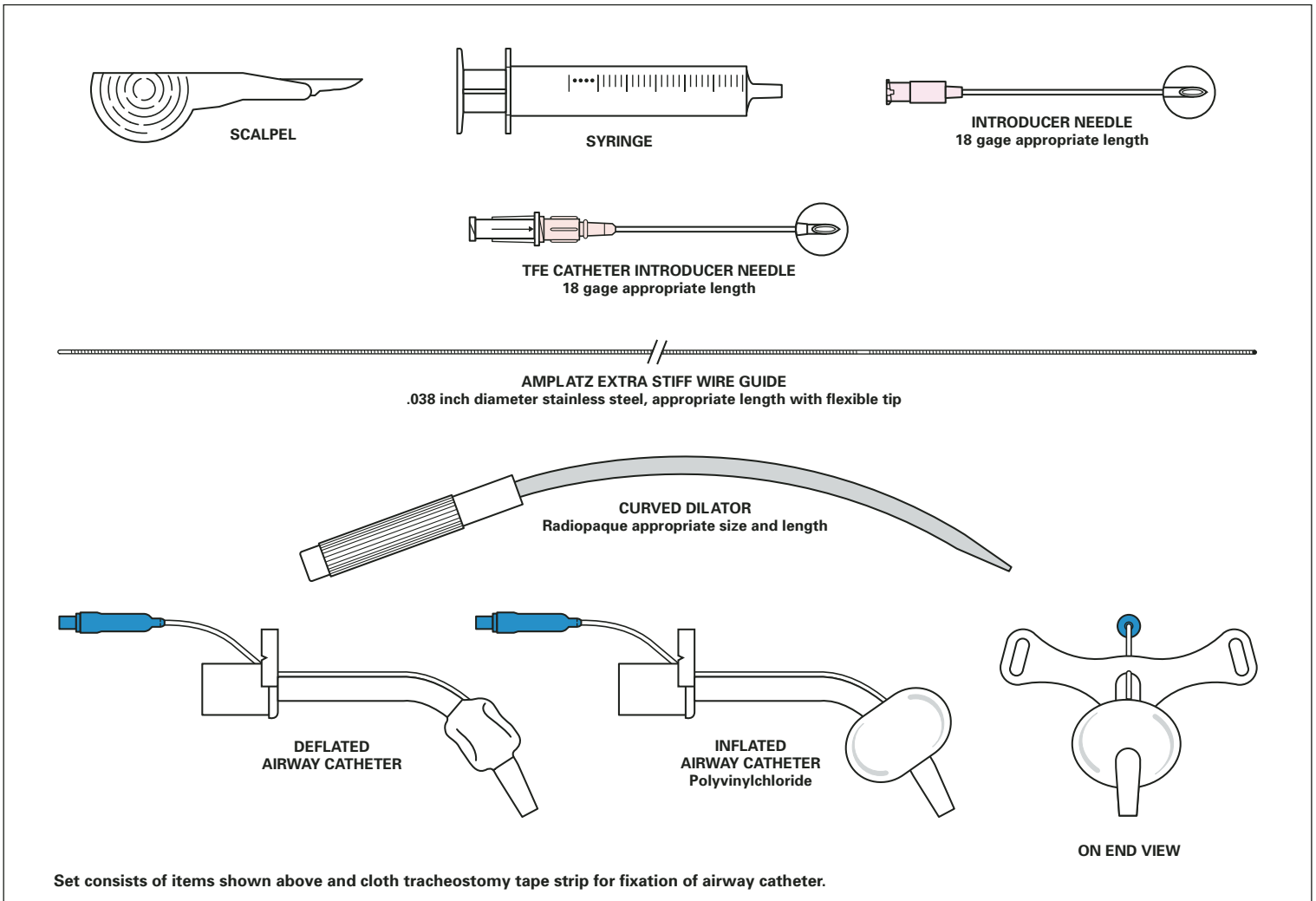

SET ORDER NUMBER	AIRWAY CATHETER		Global Product Number	Remarks
	Inner Diameter	Length		
C-TCCSB-500	5 mm	9 cm	G29432	
C-TCCSB-500-SPOPS	5 mm	9 cm	G26924	Custom packaged in slip peel-pouch design for easy transportation. Dilator is pre-inserted into airway catheter. Includes two standard, stainless steel 18 gage introducer needles, 5 cm and 7 cm long.

COOK®

Interventional Critical Care Products

www.cookcriticalcare.com

Distributed by:

COOK MEDICAL INCORPORATED

P.O. Box 4195, Bloomington, IN 47402-4195 U.S.A.

Phone: 812 339-2235, Toll Free: 800 457-4500, Toll Free Fax: 800 554-8335

COOK (CANADA) INC.

111 Sandiford Drive, Stouffville, Ontario, L4A 7X5 CANADA

Phone: 905 640-7110, Toll Free: 800 668-0300

WILLIAM A. COOK AUSTRALIA PTY. LTD.

Brisbane Technology Park, 12 Electronics Street, Eight Mile Plains

Brisbane, QLD 4113 AUSTRALIA, Phone: +61 7 38 41 11 88

WILLIAM COOK EUROPE ApS

Sandet 6, DK-4632, Bjaeverskov, DENMARK, Phone: +45 56 86 86 86