

MELKER

UNIVERSAL EMERGENCY CRICOTHYROTOMY CATHETER SET



ONE SET

TWO CHOICES

COOK®



MELKER

UNIVERSAL EMERGENCY CRICOTHYROTOMY CATHETER SET

The Melker Universal Cricothyrotomy Catheter Set features devices for both surgical and Seldinger techniques that place freedom of choice and control into the clinician's hands.

Designed to establish emergency airway access via the cricothyroid membrane for "cannot intubate, cannot ventilate" situations.



SURGICAL TECHNIQUE

This method incorporates a soft, blunted dilator with the Melker Cuffed Cricothyrotomy catheter, following the creation of the tracheal stoma.

SAFETY SCALPEL

Creates a precise incision for creation of stoma. Safety feature may help protect physician from accidental sticks.

TRACHEAL HOOK

Grasps and elevates the cricothyroid cartilage, allowing more control during airway placement.

TROUSSEAU DILATOR

Dilates by spreading stoma.

CURVED BLUNT DILATOR

Blunt-tipped dilator prevents snagging at stoma site and post-tracheal wall damage.



MELKER CUFFED EMERGENCY CRICOTHYROTOMY CATHETER



PROVEN CONTROL

Sturdy, yet flexible enough for rapid insertion, the dilator handle helps prevent slippage during catheter placement.




DEFLATED CUFF

When deflated, the cuff features a low-profile design for a smooth, atraumatic entry.



INFLATED CUFF

Low-pressure, high-volume cuff creates an airway seal to control possible aspiration and to protect the airway.

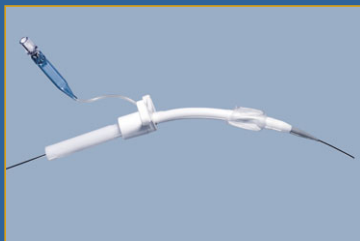

SELDINGER TECHNIQUE

Utilizes an entry needle, wire guide and dilator for a percutaneous insertion of the airway catheter.

TAPERED DILATOR

Creates gradual transition to airway catheter, enabling a smooth, easy and atraumatic insertion.

AMPLATZ WIRE GUIDE .038 INCH

Extra stiff with floppy tip.

TFE CATHETER NEEDLE

18 Gage

INTRODUCER NEEDLE

18 Gage

12 CC SYRINGE

SCALPEL

theairway
course

www.theairwaysite.com • 1-866-9AIRWAY (1-866-924-7929)

COOK®

Interventional Critical Care Products

www.cookcriticalcare.com

Distributed by:

COOK MEDICAL INCORPORATED

P.O. Box 4195, Bloomington, IN 47402-4195 U.S.A.

Phone: 812 339-2235, Toll Free: 800 457-4500, Toll Free Fax: 800 554-8335

COOK (CANADA) INC.

111 Sandford Drive, Stouffville, Ontario, L4A 7X5 CANADA

Phone: 905 640-7110, Toll Free: 800 668-0300

WILLIAM A. COOK AUSTRALIA PTY. LTD.

Brisbane Technology Park, 12 Electronics Street, Eight Mile Plains

Brisbane, QLD 4113 AUSTRALIA, Phone: +61 7 38 41 11 88

WILLIAM COOK EUROPE ApS

Sandet 6, DK-4632, Bjaeverskov, DENMARK, Phone: +45 56 86 86 86

© COOK INCORPORATED 2004

C-TCCUNI404