



20 Veterans Drive, Chicopee, MA 01022 | Tel.800.669.6900 | Fax.800.757.8585 | info@ironduck.com | www.ironduck.com

**CAPACITY:**

Pediatric Patients Aged 0-5 Years

**DIMENSIONS:**

48" x 12" x 1.75"  
Weight: 7.5lbs.

**X-RAY TRANSLUCENT:**

100% X-Ray Translucent,  
MRI and CT Scan compatible

**DESIGN:**

- Patented Dual Plane Headrop System
- Rotomolded one piece seamless HDPE shell over non-toxic Eco-Friendly polyurethane foam\*(1)

**OPTIONS:**


- P.P.E. Anti-Bac Additive
- Molded-In Graphics or Custom Screen Printing\*(2)\*(3)

**WARRANTY:**

Lifetime Warranty

**STANDARD COLORS:**

Yellow, Teal

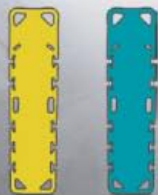

\*Pedi~Air~Align Shown with Straps.

**Pedi~Air~Align**

**Product#: 35840 (Complete)**

**Product#: 35841 (Board Only)**

The Pedi-Air-Align's patented dual plane head drop system allows for a patient to be completely immobilized while managing their airway and establishing neutral spinal alignment simultaneously. The board accommodates the larger occipital area of a child's head, allowing for neutral spinal alignment and an open, maintainable airway. The Pedi-Air-Align is the only pediatric immobilization board that is rotomolded from high-density polyethylene (HDPD) creating a one piece, seamless product with a subtle surface texture that prevents slippage but is completely impervious, and easy to clean and disinfect.

The Pedi-Air-Align Complete System includes: 1 board, 1 shoulder restraint, 1 chest strap, 4 body straps, 2 wrist restraints, 6 disposable head blocks, 6 disposable head block straps, and 1 carry case. Board size is 48" x 12" x 1.75". Lifetime guarantee.



Patented Dual Plane Headrop System.

\*(Figure 1)



**Now Available**

**P.P.E.**

- P - Protects Patients
- P - Protects EMS Professionals
- E - Extends the Life of the Board

Anti-Staph, Anti-E Coli, Anti-Pneumonia... Anti-Bac



\*(1) Kyoto & Montreal Compliant, \*(2)Free on the Purchase of 25 or more Boards \*(3)No Screen Printing on Black Backboards.

© Copyright 2009. Iron Duck, Built for Life, and Iron Duck logo are registered trademarks of Fleming Industries Inc. All Rights Reserved.